

Dry Box



Portable drying box for hearing aids

Optimal care for hearing aids as moisture and ear wax can impair function of hearings aids.

Functions:

- Special gentle drying cycle : 6 min UV-C → 3 hours 48 min drying → 6 min UV-C
- Sterilization by UV-C light (253,7 nm)
- Eliminates 99,9% of germs, approved by accredited laboratory (Bureau Veritas)
- Drying by gentle air circulation, no impact to material of your hearing heads
- Protects your hearing aids by preventing moisture
- Extends life time of your hearing aids
- Easy single-button operation
- Auto turn-off
- LED indicator for cycle
- Carrying pouch for easy transportation


Technical specifications:

- Power adapter, AC input 100-240V 50/60 Hz, DC output 5V 1A
- UV-C light (253,7 nm)
- Compartment temperature : 50° C +/- 5° C

Dimensions (width x depth x height):

- 125 x 90 x 35 mm

Content

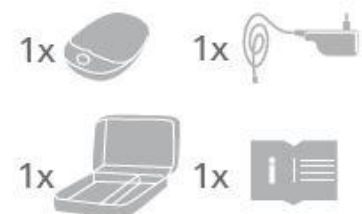

Warranty: 24 months

EAN No.: 4250711 907846

Order No.: 907846

Export Carton: 4 pcs with EAN Barcode

Picture similar / technical data / features and design are subject to permanent technical development and could be changed without prior notice – errors accepted.